

Efficacy and Safety of Novel, Disposable Endoscopic Scissors in Patients with Roux-en-Y Gastric Bypass: A Single-Center Feasibility Study

Ahmad Najdat Bazarbashi, MD¹, Thomas McCarty, MD¹, Lolwa Al-Obaid, MD², Matthew Skinner MD³, Christopher C. Thompson MD, MHES¹

¹ Division of Gastroenterology, Hepatology and Endoscopy. Brigham and Women's Hospital, Boston, MA

² Department of Internal Medicine, Lahey Hospital and Medical Center, Burlington, MA

³ Division of Gastroenterology, New York University, New York, NY

Running Title: Efficacy and Safety of Novel, Disposable Endoscopic Scissors

Correspondence to:

Christopher C. Thompson, MD, MSc
Division of Gastroenterology, Hepatology and Endoscopy
Brigham and Women's Hospital
75 Francis Street
Boston, MA 02115, USA
Phone: (617) 732-6389
Fax: (617) 264-6342
Email: cthompson@hms.harvard.edu

Word Count: 1850 (excluding abstract, references, and tables)

Figures: 3

Tables: 1

References: 21

Potential Conflicts of Interest:

Ahmad Najdat Bazarbashi has no conflicts to disclose.

Thomas R. McCarty has no conflicts to disclose.

Lolwa Al-Obaid has no conflicts to disclose.

Matthew Skinner has no conflicts to disclose.

Christopher C. Thompson is a consultant for Boston Scientific, Olympus America, and Apollo Endosurgery

ABSTRACT**Background and Aims:**

Retained suture material from primary Roux-en-Y gastric bypass or revisions may result in local inflammation, ulcer formation, and abdominal pain. The primary aim of this study was to evaluate the feasibility, efficacy, and safety of novel, disposable endoscopic scissors for suture removal.

Methods:

This was a single-center, retrospective analysis of prospectively collected data from December 2016 to January 2019. Patients with history of Roux-en-Y gastric bypass and upper endoscopy utilizing this novel, disposable endoscopic scissors device were reviewed. Measured outcomes included endoscopy indication, technical success (as determined by ability to achieve adequate cut and suture removal), improvement in abdominal pain if present prior to procedure, and adverse events.

Results:

Sixty-two patients were included in this analysis. Mean age was 54.69 ± 9.40 years. Eighty-eight percent of patients were female. Roux-en-Y gastric bypass occurred on average 142.43 ± 62.33 months prior to suture removal. Primary indications for endoscopy were evaluation of gastric pouch and gastrojejunal anastomosis for weight regain (37.10%) and abdominal pain (30.65%). Overall, technical success of these novel endoscopic scissors was 100%. Symptom improvement post-suture removal occurred in approximately half of patients (48.65%). Post-procedure bleeding was self-limited and occurred in 2 patients (3.23%). Two patients developed transient post-procedure abdominal pain.

Conclusion:

Novel, disposable endoscopic scissors appear to be highly effective and safe for removal of suture material with high technical success and minimal adverse events.

Keywords:

Obesity; Roux-en-Y gastric bypass; Endoscopy; Endoscopic Scissors

INTRODUCTION

The prevalence of obesity has nearly tripled over the past four decades, approaching a prevalence among United States adults of approximately 40%.^{1,2} Although bariatric surgery has reached a plateau over the last 2 decades, the mainstay of weight loss surgery revolves around sleeve gastrectomy and the more traditionally utilized Roux-en-Y gastric bypass procedure.³⁻⁶ Despite its efficacy in achieving weight loss, Roux-en-Y gastric bypass can be associated with a host of potential short- and long-term complications.⁷⁻⁹ While early complications may include leak, ileus, and post-operative hemorrhage, late complications can involve weight regain, anastomotic ulcers or strictures, dumping syndrome, gastro-gastric fistula, hernias, bowel obstruction, nutritional deficiencies, and abdominal pain.⁹ Although evaluation of abdominal pain in patients with Roux-en-Y gastric bypass can be challenging, often necessitating esophgogastroduodenoscopy, the cause in many cases may not be easily elucidated.^{10,11}

While the erosion of non-absorbable foreign material such as silk into the gastric pouch after bariatric surgery (usually seen in classic two-layered gastrojejunostomy) is commonly regarded as a normal or unremarkable post-operative finding, retained sutures in patients with Roux-en-Y gastric bypass have recently emerged as a potential cause for unexplained abdominal pain.^{12,13} Furthermore, additional data suggests that exposed suture may result in marginal ulcer formation with improvement illustrated post-suture removal.¹³ Based upon previous literature, it is now believed that retained suture material from primary Roux-en-Y gastric bypass or revisions may result in local inflammation, ulcer formation, abdominal pain, and also promote bezoar formation.¹³⁻¹⁵

With retained suture identified as possible sources of unexplained abdominal pain, there remains a paucity of data regarding optimal methods for suture removal through endoscopic means. Although conventional endoscopic scissors and loop cutters are reusable, they have been associated with varying technical success and may incur large costs given the need for sterilization and treatment. Recently, a novel, disposable endoscopic scissor (Ensizor, Slater Endoscopy, Miami Lakes, Florida) has been developed. Despite increasing use for various indications, limited data exists to assess this device's effectiveness, safety profile, and impact on patient symptoms. The primary aim of this study was to evaluate the feasibility, efficacy, and safety of a novel, disposable endoscopic scissor device to assist in the removal of retained sutures in patients with Roux-en-Y gastric bypass.

Methods:

Study Design and Patient Selection

This was a single-center, retrospective analysis of prospectively collected data from December 2016 to January 2019. All procedures were performed in an academic, university endoscopy unit. Data on endoscopic scissor use was obtained from multiple gastroenterologists with varying years of experience. This study was approved by the university Institutional Review Board. Eligible patients for inclusion were those patients with a history of prior Roux-en-Y gastric bypass surgery with or without previous endoscopic revision. All included patients in this study underwent use of this novel, disposable endoscopic scissor device. This data was obtained using endoscopic billing codes.

Endoscopy Protocol

Included patients underwent procedural sedation with fentanyl and midazolam. In rare instances, and when deemed appropriate by anesthesia, patients underwent deep sedation using propofol. All patients were consented for standard esophagogastroduodenoscopy with the potential for endoscopic suture removal if retained sutures were encountered. A standard upper endoscope (Olympus GIF-HQ190, Olympus America, Center Valley, PA) was used for all procedures. The gastric pouch was examined for evidence of ulceration or fistula. The gastrojejunal anastomosis was also examined for ulcers, erosions, or retained suture material.

Endoscopic Scissors and Suture Removal Technique

Once suture material was identified, the suture material was cut using the flexible endoscopic scissors device (Ensizor, Slater Endoscopy, Miami Lakes, Florida) – **Figure 1, Figure 2**. These disposable endoscopic scissors allow for the ability to cut through tissue, sutures, and many other non-metal materials using lateral forces with reinforced blades to prevent malleation, deviation, or bending. With these scissors, each blade is two components assembled together which may provide more consistent results for suture cutting. In certain instances, grasping forceps were used to pull suture material to allow more exposure at the site of mucosa attachment; thus, allowing more suture to be cut with the endoscopic scissor device. Suture material was detached from the mucosa post-incision, cut, and then retrieved with the assistance of endoscopic grasping forceps.

Outcome Measures

Primary measured outcomes of this study included technical success as determined by ability to achieve adequate cut and suture removal and adverse events. Secondary outcomes included date of Roux-en-Y gastric bypass surgery, indication for esophagogastroduodenoscopy, record of specific gastroenterology proceduralist, and assessment for improvement in abdominal pain if present prior to procedure at follow-up.

Statistical Analysis

Statistical analysis was performed using the Stata 15.0 software package (Stata Corp LP, College Station, TX). Continuous variables were expressed as means, categorical data as numbers and percentages.

Results:

Patient Demographics and Procedure Characteristics

A total of 62 patients were enrolled in this study. The mean age of included patients was 54.69 ± 9.40 years with an overwhelming percentage of females represented (n=59; 88.71%). A majority of endoscopic procedures involving suture removal were among outpatients (77.42%). The remaining cases occurred among patients admitted to the hospital. Bariatric surgery with Roux-en-Y gastric bypass procedure occurred on average 142.43 ± 62.33 months prior to suture removal. A summary table of patient characteristics is highlighted in **Table 1**.

A total of 8 different attending gastroenterologists utilized the flexible endoscopic scissors; however, a majority of cases were performed by one provider (n=50; 80.65%). The most common primary indications for upper endoscopy were for evaluation of gastric

pouch and gastrojejunal anastomosis for weight regain (n=23; 37.10%), abdominal pain (n=19; 30.65%), and nausea or vomiting (n=9; 14.52%). Other indications included: gastrojejunal anastomotic stenosis or stent exchange (n=4; 6.45%), dysphagia (n=3; 4.84%), follow-up of marginal ulceration (n=3; 4.84%), and percutaneous endoscopic gastrostomy exchange (n=1; 1.61%). Some patients had more than one indication thought the primary indication was derived from direct chart review.

Efficacy and Safety of Endoscopic Scissors Device

Surgical suture material was visualized in all included procedures with endoscopic cinches and tissue anchors noted in 4 cases. Suture associated ulceration were visualized in 29.03% of cases. Overall, endoscopic scissors were able to successfully cut the suture material in 100% of patients. Technical success of the device was 100%. Removal of cut suture material was easily accomplished with the aid of rat-tooth forceps when required.

Symptom improvement post-suture removal occurred in approximately half of patients (18 of 37 patients; 48.65%). Among the 22 patients with abdominal pain prior to procedure, 13 patients (59.09%) had resolution after suture removal. Post-procedure bleeding occurred in only 2 patients (3.23%). These bleeds were self-limited and did not require further interventions. Two patients developed transient post-procedure abdominal pain with nausea and vomiting requiring admission to the hospital. Both patients improved and were discharged few days later.

DISCUSSION

This single-center, retrospective study of a novel, disposable flexible endoscopic scissors device demonstrated an impressive technical success rate of 100%, ensuring safe and effective removal of suture material via endoscopic means with low rates of adverse events. There were no cases of unsuccessful endoscopic scissor use. Given symptomatic improvement post-suture removal, this study also validates previous literature that foreign bodies eroding into the gastric pouch may cause chronic abdominal pain in patients who have undergone Roux-en-Y gastric bypass.¹³ To these authors knowledge, this study is the first of its kind in evaluating the feasibility of these novel endoscopic scissors in patients with Roux-en-Y gastric bypass.

Retained sutures have been associated with marginal ulcerations or erosions when present at the gastrojejunal anastomosis and can cause unexplained abdominal pain in some patients. This concept of eroding foreign materials and exposed suture material causing local inflammation, ulcer formation, abdominal pain, and also bezoar formation has been demonstrated in previous literature to date.¹³⁻¹⁵ Ryou et al evaluated 21 patients with chronic abdominal pain for which the cause was believed due to exposed sutures, of which 5 patients were noted to have marginal ulcers.¹³ All patients had suture material removed with a variety of instruments including grasping forceps, loop cutters, and rat tooth forceps. Seventy-one percent of patients in this previous study reported immediate improvement of symptoms at follow-up.

Current approaches to retained sutures or foreign body material entail conventional, through-the-scope grasping forceps. This traditionally involves grasping the end of the suture and applying force by pulling the forceps back until the suture separates from the mucosa. This may often result in bleeding and potential worsening of pain.

Although this has been postulated to potentially aid in confirming the etiology through endoscopic tugging to assess reproducibility of symptoms, and identify potential individuals that would derive clinical benefit, this has not been validated in previous clinical studies.¹³ Alternatively, endoscopic through-the-scope loop cutters (Olympus America, Center Valley, PA) can be utilized. However, it is important to note that this device is not disposable, can be costly, and may carry risk of infections if not cleaned and disinfected per standard distributor recommendations.

While multiple prior studies have assessed various through-the-scope scissor-like devices, many of these have been premised on tissue resection, such as those utilized to assist in endoscopic submucosal resection using an SB-knife (SB Knife Jr; Olympus America, Center Valley, Pa, USA) or for the sole purpose of removing eroded prosthesis like vertical gastric bands.¹⁶⁻²¹ To our knowledge, this is the largest study to date assessing use of endoscopic scissors for successful resection and removal of retained sutures in patients with bariatric surgery. We highlight that the technical success of 100% was achieved despite the device being utilized by 8 different endoscopists with varying years of expertise; thus confirming its feasibility and ease of use. No appreciable learning curve was needed. The high success rate associated with these scissors may be attributed to a unique design. When compared to conventional scissors that rely upon longitudinal forces which may diminish cutting quality, this device creates a lateral force to affect cutting. This may effectively allow for improved efficacy as related to the redirected 90 degree force vector though head to head comparator studies are lacking.

We acknowledge a few limitations to this study including its single-center and retrospective design. Although this study adds to available literature that suture removal

can assist with pain and ulcer resolution as symptom improvement was noted in almost half of the patients, many of these patients were also placed on proton pump inhibitors and sucralfate which may have confounded these results. Additionally, there remains a lack of comparator studies to assess these novel, disposable scissors compared to other available products on the market. Finally, while we believe the disposable scissors do incur less costs, a cost-utility analysis would be required to compare this device to other standard methods for suture and foreign body removal.

Overall, this study reveals that novel, disposable endoscopic scissors appear to be a highly effective and safe device. The observed high technical success rate and minimal adverse events suggests these scissors may aid in removal of foreign material within the gastrointestinal tract. This notion should have special relevance to an increasing number of gastroenterologists encountering bariatric surgery patients in clinical practice.¹³ Larger studies are needed to reproduce this data and compare clinical outcomes and costs to non-disposable devices. Additionally, further studies should assess the feasibility of this device for removal of other foreign material.

REFERENCES

1. WHO. Obesity and Overweight. <https://www.who.int/news-room/fact-sheets/detail/obesity-and-overweight>. Accessed 2 April 2019.
2. Hales CM, Fryar CD, Carroll MD, Freedman DS, Ogden CL. Trends in Obesity and Severe Obesity Prevalence in US Youth and Adults by Sex and Age, 2007-2008 to 2015-2016. *JAMA* 2018;319:1723-5.
3. Khan S, Rock K, Baskara A, Qu W, Nazzal M, Ortiz J. Trends in bariatric surgery from 2008 to 2012. *Am J Surg* 2016;211:1041-6.
4. Khorgami Z, Shoar S, Andalib A, Aminian A, Brethauer SA, Schauer PR. Trends in utilization of bariatric surgery, 2010-2014: sleeve gastrectomy dominates. *Surg Obes Relat Dis* 2017;13:774-8.
5. McCarty TR, Echouffo-Tcheugui JB, Lange A, Haque L, Njei B. Impact of bariatric surgery on outcomes of patients with nonalcoholic fatty liver disease: a nationwide inpatient sample analysis, 2004-2012. *Surg Obes Relat Dis* 2018;14:74-80.
6. Jirapinyo P, Thompson CC. Endoscopic Bariatric and Metabolic Therapies: Surgical Analogues and Mechanisms of Action. *Clin Gastroenterol Hepatol* 2017;15:619-30.
7. Encinosa WE, Bernard DM, Du D, Steiner CA. Recent improvements in bariatric surgery outcomes. *Med Care* 2009;47:531-5.
8. Longitudinal Assessment of Bariatric Surgery C, Flum DR, Belle SH, et al. Perioperative safety in the longitudinal assessment of bariatric surgery. *N Engl J Med* 2009;361:445-54.
9. Ma IT, Madura JA, 2nd. Gastrointestinal Complications After Bariatric Surgery. *Gastroenterol Hepatol (N Y)* 2015;11:526-35.
10. Schulman AR, Thompson CC. Abdominal Pain in the Roux-en-Y Gastric Bypass Patient. *Am J Gastroenterol* 2018;113:161-6.
11. Dallal RM, Bailey LA. Ulcer disease after gastric bypass surgery. *Surg Obes Relat Dis* 2006;2:455-9.
12. Huang CS, Forse RA, Jacobson BC, Farraye FA. Endoscopic findings and their clinical correlations in patients with symptoms after gastric bypass surgery. *Gastrointest Endosc* 2003;58:859-66.
13. Ryou M, Mogabgab O, Lautz DB, Thompson CC. Endoscopic foreign body removal for treatment of chronic abdominal pain in patients after Roux-en-Y gastric bypass. *Surg Obes Relat Dis* 2010;6:526-31.
14. Frezza EE, Herbert H, Ford R, Wachtel MS. Endoscopic suture removal at gastrojejunal anastomosis after Roux-en-Y gastric bypass to prevent marginal ulceration. *Surg Obes Relat Dis* 2007;3:619-22.
15. Pratt JS, Van Noord M, Christison-Lagay E. The tethered bezoar as a delayed complication of laparoscopic Roux-en-Y gastric bypass: a case report. *J Gastrointest Surg* 2007;11:690-2.
16. Ge PS, Thompson CC, Aihara H. Endoscopic submucosal dissection of a large cecal polyp using a scissor-type knife: implications for training in ESD. *VideoGIE* 2018;3:313-5.

17. Oka S, Tanaka S, Takata S, Kanao H, Chayama K. Usefulness and safety of SB knife Jr in endoscopic submucosal dissection for colorectal tumors. *Dig Endosc* 2012;24 Suppl 1:90-5.
18. Yamashina T, Takeuchi Y, Nagai K, et al. Scissor-type knife significantly improves self-completion rate of colorectal endoscopic submucosal dissection: Single-center prospective randomized trial. *Dig Endosc* 2017;29:322-9.
19. Evans JA, Williams NN, Chan EP, Kochman ML. Endoscopic removal of eroded bands in vertical banded gastroplasty: a novel use of endoscopic scissors (with video). *Gastrointest Endosc* 2006;64:801-4.
20. Shehab H, Gawdat K. Endoscopic Management of Eroded Bands Following Banded-Gastric Bypass (with Video). *Obes Surg* 2017;27:1804-8.
21. Spann MD, Aher CV, English WJ, Williams DB. Endoscopic management of erosion after banded bariatric procedures. *Surg Obes Relat Dis* 2017;13:1875-9.

TABLE**Table 1.** Summary Study Characteristics of Roux-en-Y Patients with Suture Removal

Number of Patients	62
Mean Age (years)	54.69 ±9.40
Number of Females (%)	59 (88.71%)
Time Since Roux-en-Y (months)	142.43 ±62.33
Primary Indications for EGD	
Weight Regain	23 (37.10%)
Abdominal Pain	19 (30.65%)
Nausea and Vomiting	9 (14.52%)
Gastrojejunal Stenosis	4 (6.45%)
Dysphagia	3 (4.84%)
Follow-up of Marginal Ulcer	3 (4.84%)
PEG exchange	1 (1.61%)

EGD: esophgogastroduodenoscopy

PEG: percutaneous endoscopic gastrostomy exchange

FIGURES

Figure 1. Gross appearance of novel, disposable endoscopic scissors.



Figure 2. Unique blade design of novel, disposable endoscopic scissors.

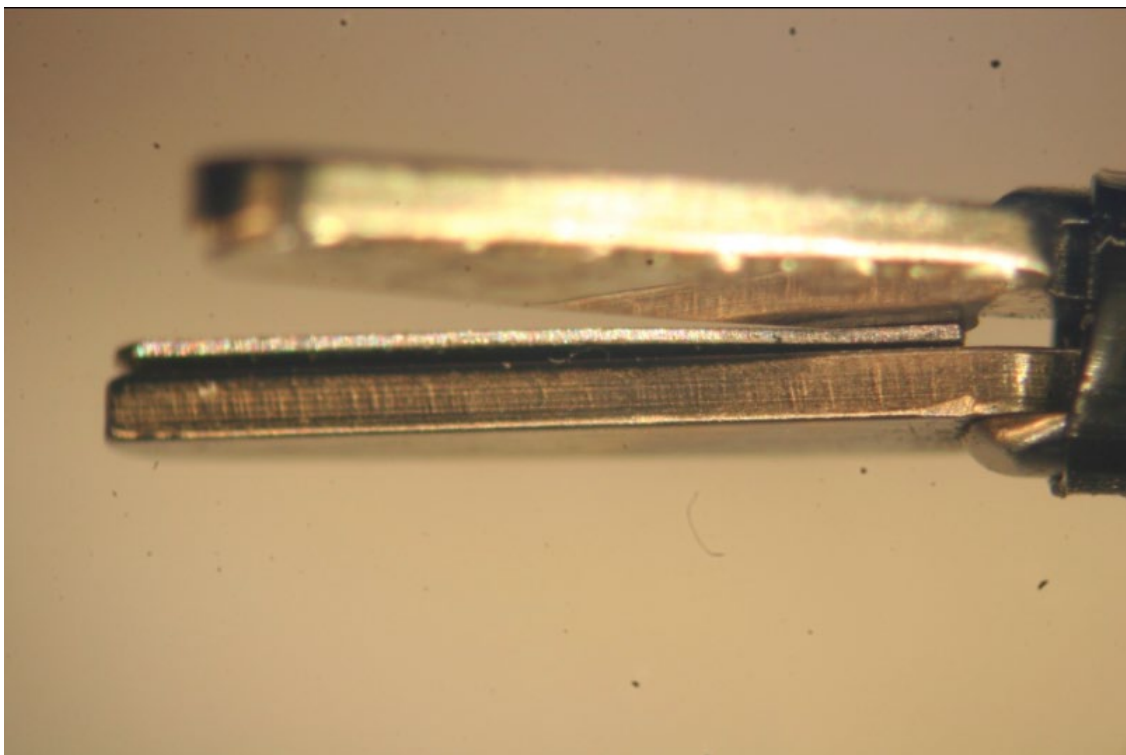


Figure 3. Esophgogastroduodenoscopy of patient with Roux-en-Y gastric bypass visualizing suture associated ulceration and successful removal with endoscopic scissors.

